


LikoLight

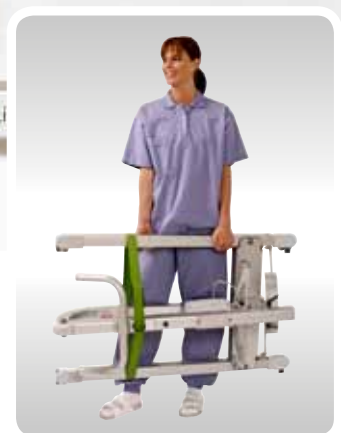

The LikoLight is developed for those in need of a light mobile lift that folds away easily. LikoLight works much like a normal electric mobile lift, but the lift is portable and can be conveniently stored. In fact, LikoLight can be quickly disassembled into three pieces to facilitate transport and storage.

Each detail on the LikoLight is carefully designed to provide maximum stability but minimum weight. The total weight is only 24 kg (53 lbs.) and no part weighs more than 10 kg (22 lbs).

Practical LikoLight Travel Bags are available as an accessory. The bags are manufactured in a padded, durable light-weight material, and are equipped with wheels for easy transport.

The versatile lightweight lift

Despite its size and weight, LikoLight is ideal for all common lifting situations; for example, lifts from bed to wheelchair, to and from the toilet and bathtub, or for lifting someone who has fallen to the floor.

LikoLight has electric raising and lowering and mechanical base-width adjustment. A patented easy-to-use mechanical emergency lowering is standard.

LikoLight Travel Bags are practical carrier bags made of a padded, durable and soft lightweight material. Perfect for travels.

Important features:

- Folds away easily.
- Easily disassembled into three pieces.
No part weighs more than 10 kg (22 lbs.).
- Made of aluminium; the right choice for environment, quality and design.
The choice of material also makes the lift easy to roll and easy to keep clean.
- Mechanical emergency lowering.
- Stable construction.
- Ideal for travelling.
- Practical tabletop charger available as an accessory.



Accessory: LikoLight Travel Bags.



Easy to bring on the trip.



Manages all common lifting situations, including lifts from the floor.

■ TECHNICAL DATA

Lifting capacity:	140 kg (308 lbs., 22 st.)
Lifting speed:	3.8 cm/s (1.5 inch/s) without load
Batteries:	2 x 12 V 2.9 Ah, valve regulated, sealed, lead accumulator (gel-type batteries)
Battery charger:	100-240 V 50-60 Hz, max. 400 mA
Material:	Aluminium
Weight:	24 kg (53 lbs.) total weight; disassembled, no part weighs more than 10 kg (22 lbs.)
Product numbers:	2030001 LikoLight 2037025 LikoLight Travel Bags
Patented.	



LikoLight is tested by an accredited testing institute and complies with the requirements of directives for medical technical Class 1 products (MDD 93/42/EEC).

LikoLight complies with the requirements according to IEC 60601-1, IEC 60601-1-2, SS-EN ISO 10535, UL-2601-1 and CAN/CSA C22.2 No.606.1.

Made in Sweden.

Liko is quality certified according to ISO 9001 and its equivalence for the medical device industry, ISO 13485.

Liko is also certified according to environmental standard ISO 14001.



www.liko.com

Manufactured by:

Liko AB
Nedre vägen 100
975 92 Luleå
Sweden
info@liko.se

Local representative in UK

Hill-Rom Limited
Clinitron House, Ashby Park
Ashby de la Zouch,
Leicester LE65 1JG
Tel: 01530 411000
Fax: 01530 411555
info@liko.co.uk

To find your local Liko representative,

please visit www.liko.com or call
+46 (0)920 47 47 00



A **Hill-Rom** Company